



Thank you for purchasing this product. Please read this manual carefully before using and keep well.

1 Safety precautions - For safety use -

Please read this Safety precautions carefully.

For safety using, please note below cautions that represents of the harm or damage caused by improper handling.
If you give or lend this product to others, note to give with this Safety precautions together.

[Display means]

Display	Means
	The wrong way to use makes people die or risk of serious injury.
	The wrong way to use be possible to make people die or risk of serious injury.
	The wrong way to use cause to damage to person or property.

[Picture display]

	danger/warning/attention
	Prohibited matter
	Attention and instruction

Before use

⚠ Danger

- People who use these following medical electronics, do not use this machine
 - Implantable medical electronic equipment
 - Medical electronic equipment life support
 - Medical electronic equipment mounted type
 - Other medical electronic equipment
- ※ They may have cause the error.

⚠ Warning

- People who applies to the following
 - Acute disease
 - Febrile disease
 - Tuberculous diseases
 - Hemophilia disease
 - Abnormal blood pressure
 - Malignant tumor
 - Infectious disease
 - Heart disease
 - Angina pectoris
 - Multiple sclerosis
 - Purpura, hemorrhage
 - Skin cancer
 - Skin abnormalities
 - Pregnancy
 - Lactation
 - Epilepsy
 - Facial neuralgia
 - Keloid diathesis
 - Atopy
 - Excessive tanning
 - Purulent inflammation of acne
 - Allergic constitution
 - Rash experience
 - Steroid hormone
 - Precision eye disease
 - Infants and young children
 - People who cannot control themselves
 - During menstruation
 - Other under treatment and abnormality in the body

Place/Power

⚠ Attention

- This product is not water-proof. Please do not use or charge in bathroom.
- Do not put it in a place where fire and high temperature & a place where direct sunlight.
- Do not use plug and code in dirty place.

⚠ Attention

- Please use the specified conditions (voltage, cable etc). Otherwise it will cause fire, electric shock.
- In addition when pull out the plug from power supply, do not pull the code. It will cause the electric contact failure.

Prepare and use note

⚠ Attention

- Please use on the face.
- Please use the method of the safety precaution.
- Please follow the usage of time.

- This product is domestic electric appliance, do not use for commercial purposes.
- If the product and the code have any damage, please stop using it.
- Do not touch the machine and the code with wet hands.

⚠ Warning

- If you feel sick, please stop using and consult to your doctor immediately.
- If the product has any problem, please stop using.

⚠ Attention

- Please follow the below rules
 - Make sure that there is no abnormality feeling after treatment.
 - Put off the all metal products and accessories.
 - Don't apply on hair removed skin.
 - ※ If it appears any irritations, please consult your doctor.
 - Don't use the product over proper time.
 - Make sure that all function working before using.

Repair-Maintenance

⚠ Attention

- Do not put the product in places where children can touch.
- Please turn off the power when you clean the product.

- Do not transform or take the product apart.


⚠ Warning

- Do not put in high temperature or direct sunlight place.
- Do not dispose this product like ordinary household waste when the expected life is over. Please dispose this product according to the local government regulation.

- Please do not use to clean the product with diluent, essential oil, gasoline, containing drugs, cleaning agents, and hot water.
- Please do not put the product around the hot machine.
- Please do not fall the product.
- Please do not wash with water.

-1-


2 Set *Accessories and packaging are subject to change

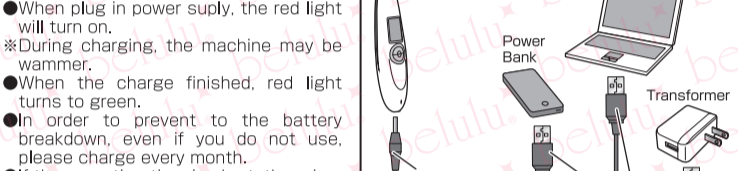

[Product Specification]
Charging Time: 4~8hours Input: 5V/0.5A
Battery type: Lithium-Ion
Dimensions: 180x45x35(mm) Weight: 100g
[AC adapter specification]
Power Input: 100~240V 50/60Hz 0.2A Power Output: 5V/1A

3 Charge

- Please charge at first.
- When plug in power supply, the red light will turn on.
- ※ During charging, the machine may be warmer.
- When the charge finished, red light turns to green.
- In order to prevent to the battery breakdown, even if you do not use, please charge every month.
- If the operating time is short, there is a possibility that the battery life is over.

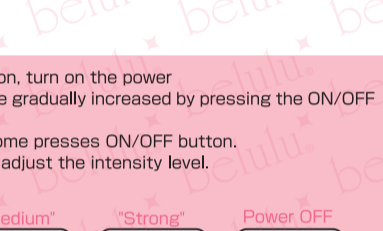


-2-

4 Operation

ON/OFF button EMS/RF


- Press the ON/OFF button, turn on the power
- The intensity level will be gradually increased by pressing the ON/OFF button each time.
- Turn off the power by some presses ON/OFF button.
- Restart the machine to adjust the intensity level.

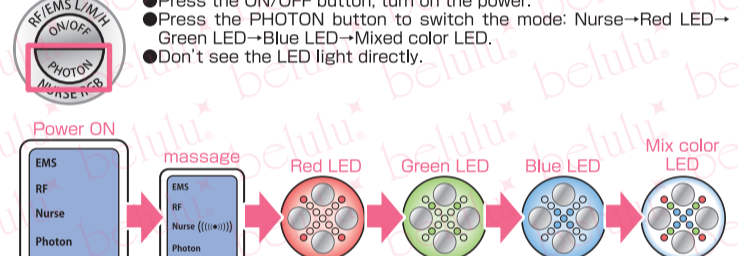

● You can use Nurse massage or LED photon at the same time.

PHOTON button Message & LED photon

- Press the ON/OFF button, turn on the power.
- Press the PHOTON button to switch the mode: Nurse → Red LED → Green LED → Blue LED → Mixed color LED.
- Don't see the LED light directly.

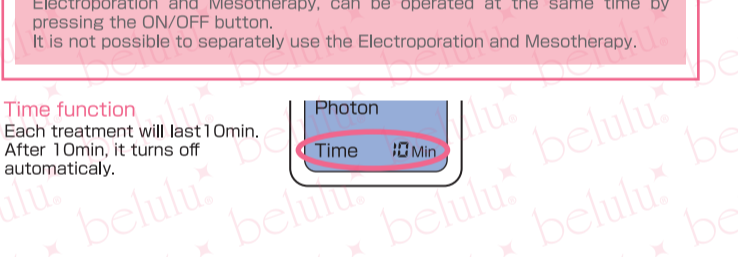

● You can use EMS/RF mode at the same time.

Combination

- EMS/RF Treatment + Massage
 - EMS/RF Treatment + LED(PHOTON)
- Electroporation and Mesotherapy, can be operated at the same time by pressing the ON/OFF button.
It is not possible to separately use the Electroporation and Mesotherapy.

Time function

Each treatment will last 10min. After 10min, it turns off automatically.



-3-

[Attention]

- Please check allergy test at first.
- If there are any problems, stop using, and consult to your doctor.
- When use the LED light, don't see the LED light directly.
- When use the ion introduction, please apply any liquid lotion.
- Please take the side metal part to improve the ion induction.
- Put the 3D head firmly on your face and softly slide when you massage.
- Please put off all accessories and other metal products.
- If you feel sick, stop using and consult to your doctor immediately.

Examples of treatment.

Facial treatment EMS/RF mode

- Apply any beauty liquid.
 - Turn on the power by pressing ON/OFF. And press again to start EMS / RF mode.
- You can add the lotion or serum to avoid the skin damage by friction.


For dry skin Electroporation/Mesotherapy+Nurse

- Apply any beauty liquid.
 - Turn on the power by pressing ON/OFF. And press again to start EMS / RF mode. Then press the PHOTON to start Nurse massage.
- You can use it on the face mask

For Dullness EMS/RF & Massage + Green LED

- Apply any beauty liquid.
 - Turn on the power by pressing ON/OFF. And press again to start EMS / RF mode. Then press the PHOTON to start Nurse massage.
 - Later, press the PHOTON button to turn on the Green LED light.
 - Start the green LED treatment
- Don't see the LED light directly. Please close your eyes.

FACE CARE Each left side 2min right side 2min, total 4 minutes

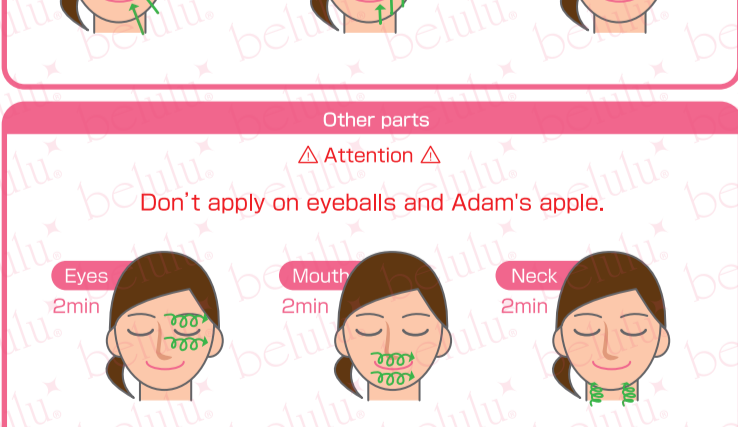

Other parts

⚠ Attention ⚠

Don't apply on eyeballs and Adam's apple.



You may feel a little electric shock when use EMS. It's not electric leakage. Please adjust the level of EMS/RF power.

[Frequency of use] EMS/RF : 3days / week LED light : everyday

5 After the treatment

- Turn off the power by a long pressing.
- Remove the all liquid on product.
- Clean 3Dhead with disinfectant alcohol.
- ※ Do not use any oil
- ※ Do not wash with water
- Do not put the product in high temperature and high moisture areas. Do not put the product where children can touch.

-4-

Q&A

Q When using EMS/RF, it doesn't vibrate.
A EMS/RF mode is just an electric stimulation and warming treatment. So it doesn't vibrate.

Q Can I use it every day?
A We recommend the usage of 3 times in a week. Or you can treat LED PHOTON everyday. Don't use too much and keep proper usage and time.

Q Should I use any beauty liquid?
A Please apply any beauty liquid to avoid skin damage.

Q What is the 3D head made of?
A The head metal parts are special steel which is expected not to be metal allergy.

Q Why does it count down from 10min?
A By pressing each button, the timer will start again. After 10 minutes, it will turn off automatically.

Q The display and LED light looks weaker.
A You need to charge it.

If you have any questions please contact us.

JAPAN
TEL : 045-227-5123
E-Mail : shop@kireido.com

CHINA
TEL : 0771-28752414
E-Mail : china@kireido.com

-5-